



## UB-328 Wrist Auto-Inflation



**One of the most technologically advanced yet portable blood pressure monitors available in the marketplace today.**


A&D Medical has applied its state-of-the-art measurement technology to monitoring blood pressure at the wrist. The Wrist Auto-Inflation unit is activated with the push of a single button. The large display simultaneously shows systolic, diastolic and pulse readings. The automatic memory keeps your last 30 blood pressure measurements close at hand.

This monitor stores the last 30 readings in memory and displays average readings. The average reading feature automatically calculates the “average” of the total readings stored in memory. The UB-328 also has a unique Irregular Heartbeat feature. This monitor provides blood pressure and pulse rate measurements even when an Irregular Heartbeat occurs. A rigid carrying case is included to protect the monitor when not in use.

- One button operation
- Simultaneous readout of systolic/diastolic pressure and pulse rate
- Large-print instruction guide with illustrations and quick start card in English, French and Spanish
- Four week blood pressure logbook included
- Rigid carrying case included

For questions about **LifeSource** monitors, call us toll-free at 1-800-726-3364 and a specially trained representative will assist you, or visit us at [www.LifeSourceOnline.com](http://www.LifeSourceOnline.com).

- Professional accuracy
- 30 memory recall
- Irregular Heartbeat feature
- Calculates and displays average readings
- Fast measurement
- Ideal for travel

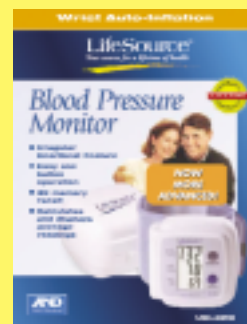

# UB-328 Wrist Auto-Inflation

The UB-328 Wrist Auto-Inflation unit straps on quickly and comfortably to one's wrist to provide an accurate reading. This compact and portable unit gives fast measurements and provides convenience to those looking for an easy way to monitor their blood pressure.

	PART NUMBER	UPC NUMBER	CUFF SIZE
UNITED STATES AND CANADA	UB-328	0-93764-60036-4	5.3" - 7.7" (13.5-19.5cm)



Each unit comes with a rigid carrying case, trilingual instruction manual, quick start card, four week blood pressure logbook and warranty card.

## Technical Specifications

<b>Dimensions</b>	3.1 (l) x 2.9 (w) x 2.3" (h) (79 x 73 x 58mm)
<b>Dimensions with Case</b>	4.0 (l) x 3.0 (w) x 3.3" (h) (102 x 76 x 83mm)
<b>Weight</b>	4.2 oz. (120 g) without batteries
<b>Display</b>	Digital: 10/7 mm character height Pressure and pulse displayed simultaneously
<b>Memory</b>	30 readings
<b>Pressurization</b>	Automatic, using micropump
<b>Deflation</b>	Constant air release-valve system
<b>Type</b>	Oscillometric
<b>Accuracy</b>	Pressure: ±3 mmHg or 2%, whichever is greater; Pulse: ±5%
<b>Pressure Range</b>	20mmHg to 280 mmHg; Pulse: 40 to 200 pulses/minute
<b>Power Source</b>	Two type AAA alkaline batteries (1.5 volt). Batteries not included.
<b>Battery Life</b>	Approx. 400 measurements
<b>Operating Environ.</b>	50° to 104°F (10° to 40°C)
<b>Storage Environ.</b>	14° to 140°F (-10° to 60°C)

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute (ANSI/AAMI SP10) for electronic sphygmomanometers. Manufactured by an ISO 9001 and ISO 13485 registered company. Engineered in Japan. Made in China.

## Packaging and Shipping Specifications

<b>Selling Unit (Box)</b>	<b>L x W x H (inches)</b>	<b>WT (lb.)</b>	<b>Individual Item Cube (cubic feet)</b>
UB-328	4 <sup>9</sup> / <sub>16</sub> x 3 <sup>7</sup> / <sub>16</sub> x 6 <sup>1</sup> / <sub>16</sub>	0.65	0.06
<b>Unit Case (10 units)</b>	<b>L x W x H (inches)</b>	<b>WT (lb.)</b>	<b>Case Cube (cubic feet)</b>
UB-328	18 x 10 x 7	7.75	0.73



Ask your A&D Sales Representative about other *LifeSource* products.

- Blood Pressure Monitors
- Heart Rate Monitors
- Nebulizers
- Digital Thermometers
- Personal Scales
- Stethoscopes

For questions about *LifeSource* monitors, call us toll-free at 1-800-726-3364 or visit our web site at [www.LifeSourceOnline.com](http://www.LifeSourceOnline.com)



1555 McCandless Drive  
Milpitas, CA 95035  
[www.LifeSourceOnline.com](http://www.LifeSourceOnline.com)

LifeSource is a trademark of A&D Medical. ©2004 A&D Medical. All rights reserved.

**LifeSource**<sup>™</sup>  
Your source for a lifetime of health



LIT:LS328 5/04