




Welch Allyn
SureTemp® Plus 692 Electronic Thermometer

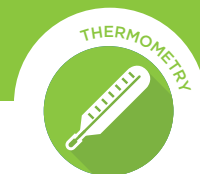

The SureTemp Plus 692 Thermometer— with the speed and accuracy you require, and the safety and security you demand.

1

Easy:

Intuitive user interface, easy-to-read LCD screen and flashing icons, helps select and identify temperature measurements and modes.

2

Safe:

Designed with interchangeable oral/axillary (blue) and rectal (red) probe wells to address patient safety and help reduce the risk of cross-contamination.

3

Fast and Accurate:

With proven accuracy, the SureTemp 692 thermometer captures quick readings—approximately 4-6 second oral, 10-13 second pediatric axillary and rectal readings and 12-15 second adult axillary readings.

4

Secure:


To top it all off, our security features make sure the unit stays where it is supposed to be. A location ID on the display and other innovative features help prevent theft.

- Easy-to-use, intuitive operation
- Interchangeable oral/axillary and rectal probe wells help reduce the risk of cross-contamination



SureTemp Plus 692—setting the standard of care

ID field for electronic labeling of thermometer

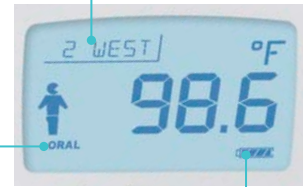

User-selectable icon and/or words for identification of thermometer mode

Battery life indicator

Interchangeable probe with push-button probe cover ejection

Removable probe well for easy cleaning

C/F button for toggling between Fahrenheit and Celsius

Prominent recall button for display of last temperature taken

Ergonomic design with rubberized grips

Convenient storage for 25 probe covers

Large, easy-to-read, backlit LCD

Mode button for selecting oral, pediatric axillary or adult axillary

Pulse timer button/monitor mode

NEW feature makes it easier to switch to Monitor Mode for measuring a patient's thermal steady state (approx. 3 to 5 minutes)

Interchangeable rectal probe and well reduce risk of cross-contamination



The **SureTemp Plus 692** with **standard wall mount** has storage for an additional 50 probe covers.



The SureTemp Plus 692—setting the standard of care even higher

SureTemp Plus 692 Product Specifications

Calibration Accuracy	±0.2 °F (±0.1 °C) (Monitor Mode)	
Temperature Predict Time	Oral: Adult Axillary (18 years and older): Pediatric Axillary (17 years and younger): Rectal:	Approximately 4-6 secs Approximately 12-15 secs Approximately 10-13 secs Approximately 10-13 secs
Unit Dimensions	8.46" x 3.18" x 2.43" (215 mm x 81 mm x 62 mm)	
Weight	12.6 oz (357 g)	
Ambient Temperature Range	50.0 °F (10 °C) to 104 °F (40.0 °C)	
Patient Temperature Range	80.0 °F (26.7 °C) to 110 °F (43.3 °C)	
Power Source	(3) 1.5Vdc AA batteries	
Battery Life	Approximately 5,000 temperatures	

Ordering Information:

THERMOMETERS

- 01692-200** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Security System with ID Location Field, Oral Probe and Oral Probe Well with 4ft cord
- 01692-201** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Security System with ID Location Field, Rectal Probe and Rectal Probe Well with 4ft cord
- 01692-300** SureTemp Plus 692 Electronic Thermometer, One per room, Wall Mount, Oral Probe and Oral Probe Well with 9ft cord (must be locked to the wall)
- 01692-301** SureTemp Plus 692 Electronic Thermometer, One Per Room, Wall Mount, Rectal Probe and Rectal Probe Well with 9ft cord (must be locked to the wall)
- 01692-700** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Oral Probe and Oral Probe Well with 9ft cord, with Rolling Stand and Holder for Spare Probe & Well

ACCESSORIES

- 21326-0000** Wall Holder
- 21330-0000** Holder for Spare Probe and Well
- 06138-000** Calibration Key Assembly, 690/692
- 02892-000** Probe and Well Kit, 4ft, Rectal
- 02892-100** Probe and Well Kit, 9ft, Rectal
- 02893-000** Probe and Well Kit, 4ft, Oral
- 02893-100** Probe and Well Kit, 9ft, Oral
- 02891-0000** Probe Well, Blue
- 02891-1000** Probe Well, Red
- 08273-000** Thermometer Stand with Holder for Spare Probe and Well
- 21333-0000** Bracket, M600 Stand Upgrade, with Holder for Spare Probe and Well
- 98025-0000** Extended Warranty M692

REPLACEMENT PROBE COVERS

- 05031-101** Case of 1,000 Probe Covers

Contact your Welch Allyn representative or visit www.welchallyn.com to learn more.



Welch Allyn, Inc.
 4341 State Street Road
 Skaneateles Falls, NY 13153 USA
 (p) 800.535.6663 (f) 315.685.3361